




# PADI4 Rabbit pAb

货号	YP-Ab-18313
同位型	IgG
应用	ELISA,WB,IHC-P
种属	Human, Mouse, Rat
免疫原	Recombinant fusion protein containing a sequence corresponding to amino acids 290-663 of human PADI4 (NP_036519.2).
来源	Rabbit
稀释	WB 1:500 - 1:2000 IHC-P 1:50 - 1:200
纯化工艺	Affinity purification
分子量	74kDa
背景	This gene is a member of a gene family which encodes enzymes responsible for the conversion of arginine residues to citrulline residues. This gene may play a role in granulocyte and macrophage development leading to inflammation and immune response.
浓度	1 mg/ml
储存	-15°C to -25°C/1 year(Do not lower than -25°C)
有关注意事项	Avoid repeated freezing and thawing!
使用建议	This product can be used in immunological reaction related experiments. For more information, please consult technical personnel.

## Products Images

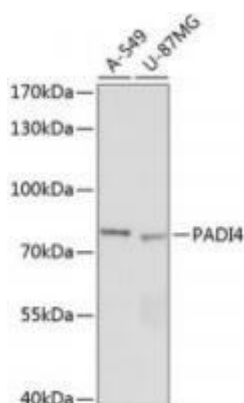


Western blot analysis of various lysates using PADI4 Rabbit pAb (A16188) at 1:1000 dilution.

Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: 25µg per lane.


Blocking bufer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 90s.

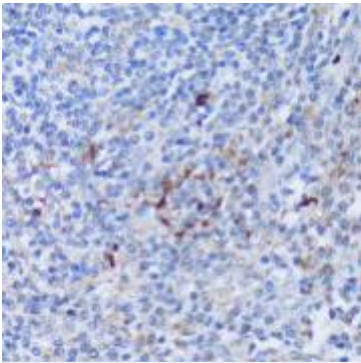


Immunohistochemistry analysis of paraffin-embedded Human tonsil using PADI4 Rabbit pAb (A16188) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

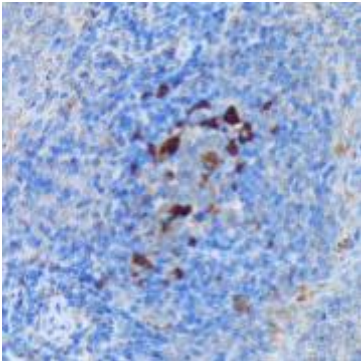

Immunohistochemistry analysis of paraffin-embedded Rat spleen using PADI4 Rabbit pAb (A16188) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

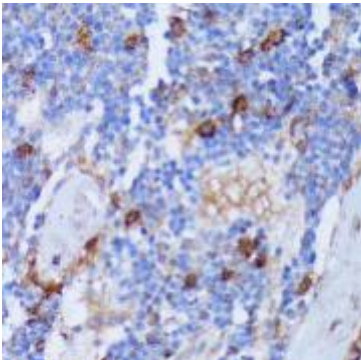

Immunohistochemistry analysis of paraffin-embedded Rat bone marrow using PADI4 Rabbit pAb (A16188) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

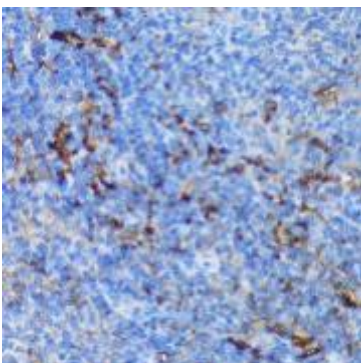

Immunohistochemistry analysis of paraffin-embedded Mouse spleen using PADI4 Rabbit pAb (A16188) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.